


Technical data and operating instructions.
For in vitro use only.

Vivapure[®] AdenoPACK[™] 20 RT

Adenovirus (Ad5) purification and concentration kit for up to 20 ml cell culture volume (E.g. 1 + 15 cm plate)

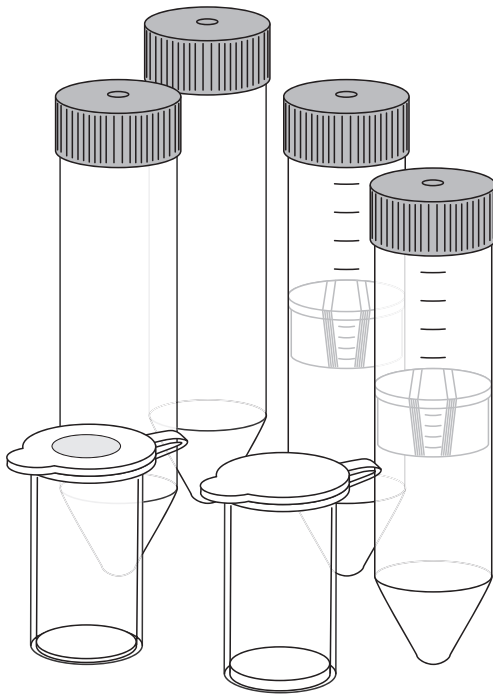

Vivapure AdenoPACK 20 RT– Introduction

Storage conditions | shelf life

The AdenoPACK kit contents should be stored at room temperature. This kit should be used within 24 months of purchase.

Introduction

This protocol describes the purification of Adenovirus (Ad5 strains) with Vivapure Q Maxi M spin columns containing an ion exchange membrane adsorber that binds adenoviral particles. Once bound, virus particles can be further purified by washing away non-specifically bound proteins, before elution within one hour.

The Vivapure Q Maxi M is designed for parallel purification and concentration of Ad5 strain adenovirus from 20 ml culture medium.

In contrast, traditional CsCl gradient centrifugation is a time consuming method, typically taking 1–2 days, which can only be performed efficiently for cell culture volumes larger than 100 ml. Furthermore, the toxicity of the media places limitations on downstream applications.

Ready-to-use filter devices, Vivapure Q Maxi M spin columns, centrifugal Vivaspin concentrators and buffers make the following purification procedure as easy as filtration.

Virus purification tests conducted in cooperation with Progen Biotechnik GmbH, Heidelberg.



Vivapure AdenoPACK 20 RT

Cat. Number	VS-AVPQ022
Number of purifications possible with AdenoPACK 20	6 × 20 ml
Vivapure Q Maxi M spin columns	6
Vivaclear Maxi 0.45 µm PES	6
Empty 50 ml tubes	6
Loading Buffer (10×)	25 ml
Washing Buffer (10×)	30 ml
Elution Buffer	20 ml
Vivaspin 20, 100 kDa MWCO	6
Instructions	1 each for Kit and Vivaspin


Materials of construction

Vivapure Q Maxi M spin column housing	Polypropylene
Vivaclear Maxi housing	Polypropylene
Vivapure Q Maxi M membrane	Stabilised RC
Buffer containers	LDPE
Purification buffers	Proprietary


Kit specifications

Sample size	20 ml of Adenovirus supernatants
Virus particles (VP) per ml	Typically up to 1×10^{11} – 10^{12}
VP/IU	50–100
Processing time	Typically 1 hour
Endotoxin level	<0.025 EU/ml

Kit contents

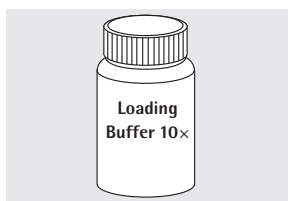


1 × 25 ml Loading Buffer

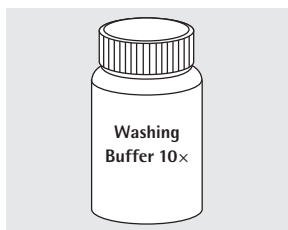


1 × 30 ml Washing Buffer




1 × 20 ml Elution Buffer

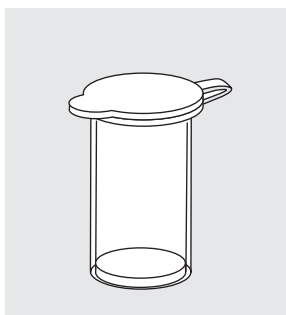


6 × Vivapure Q Maxi M spin columns

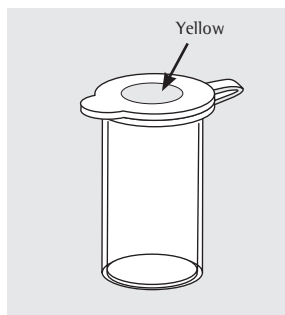

6 × Vivaclear Maxi spin columns

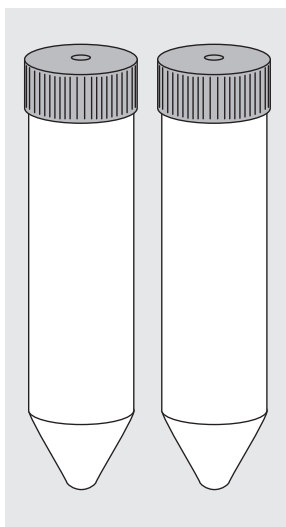

6 × 50 ml tubes with screw caps

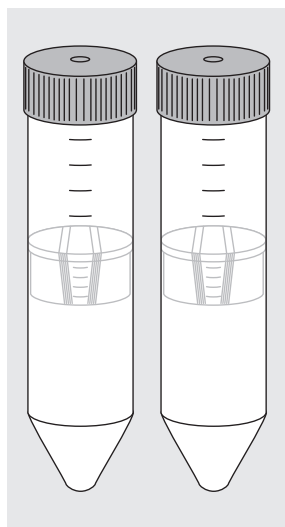

6 × Vivaspin 20 concentrators with 100 kDa membrane

Additional material required but not supplied

Centrifuge with swing-out rotors accepting 50 ml falcon tubes

Ethanol | dry ice bath or -80°C freezer


Water bath at 25°C

Sterile plastic container for sample handling

Optional – Storage Buffer: 20 mM Tris/HCl, 25 mM NaCl, 2.5% Glycerol (w/v), pH 8.0 at 22°C

1500 U Benzonase® Nuclease

Purification protocol – Overview

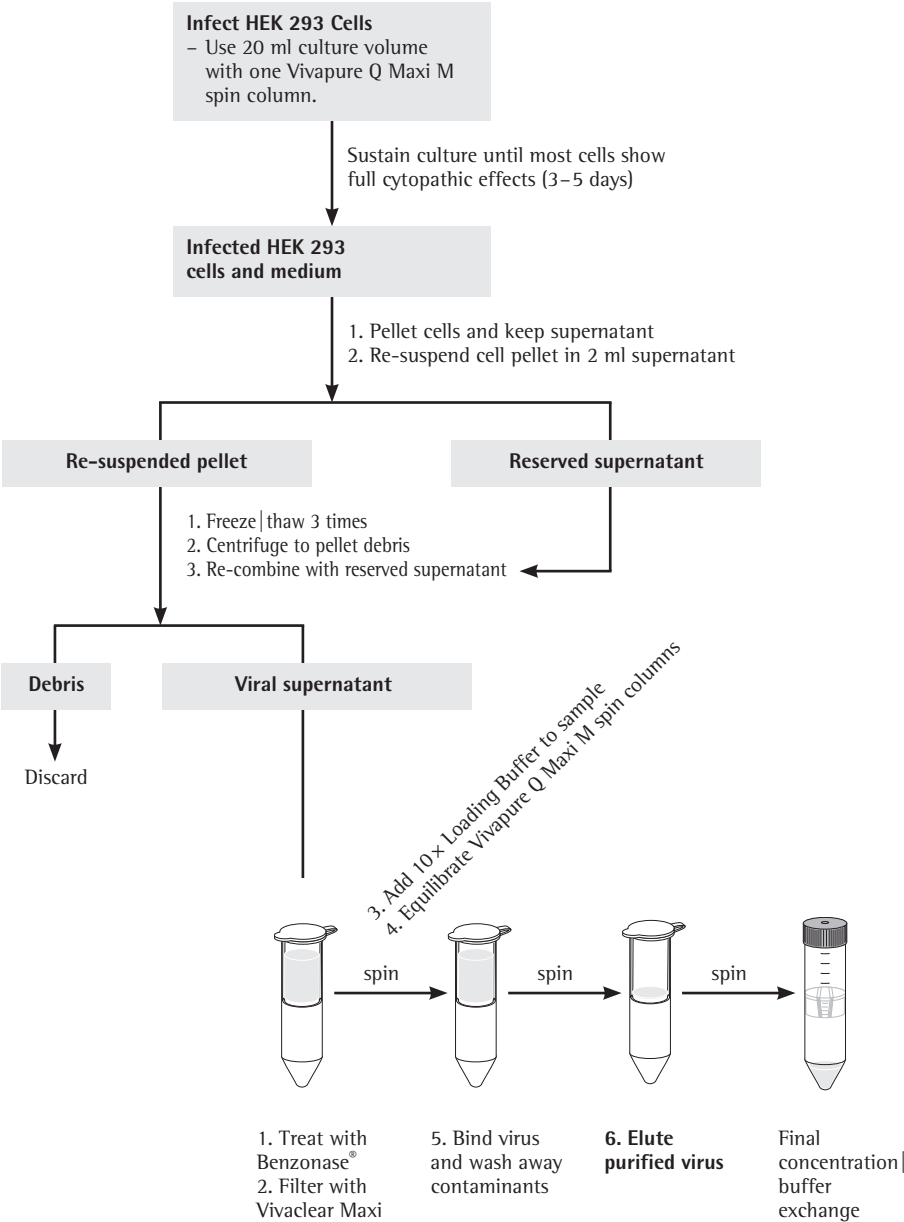

Purification protocol – Techniques

A). Sample preparation

Note:

Each Vivapure Q Maxi M spin column may be used to purify virus from up to 20 ml culture volume. This kit contains sufficient consumables for six such small preparations. Please adjust reagent volumes accordingly for smaller samples.

1. Amplify adenovirus in low-passage HEK 293 cells in up to 20 ml total culture volume (E.g. one 15 cm plate with 20 ml culture) that has been infected with an adenovirus stock at an m.o.i. of 10–20. Cultures should be grown in DMEM +10% FBS pH 7.0–7.4 at 37°C with 5% CO₂.

2. Once most of the cells show cytopathic effects (2–5 days), pool cells and medium. It may be necessary to detach adhering cells using a pipette or cell scraper.

3. Centrifuge at 3,500 × g for 15 minutes to pellet cells.

4. Decant supernatant to a sterile container and set aside.

5. Re-suspend cell pellet in 2 ml supernatant.

6. Freeze – thaw the cell suspension completely 3 times to disrupt cells alternately using a 25°C water bath and ethanol| dry-ice bath or –80°C freezer.

7. Centrifuge at 3,500 × g for 15 minutes to pellet cell debris.

8. Decant viral supernatant, re-combine it with the original supernatant and mix gently.

9. Add Benzonase® Nuclease, to a final concentration of 12.5 U/ml.

10. Mix sample and incubate for 30 min at 37°C in order to digest cellular nucleic acids.

11. Load digested supernatant on a Vivaclear Maxi and spin the device 5 minutes at 500 × g or until the whole volume has passed the membrane.

12. Collect the flow-through. Estimate the volume and slowly add 1/9 volume of 10 fold loading buffer under agitation to avoid osmotic shock in the virus particles. E.g. 2 ml to 18 ml flow-through. Accurate volume measurements of the supernatant and 10 × Loading Buffer are critical to achieving the right conditions for binding virus particles.

B). Vivapure Q Maxi M preparation

13. Dilute 10 × Washing Buffer to working concentration, e.g. for one preparation: 5 ml buffer with 45 ml deionised water, and mix well.

14. Equilibrate the Vivapure Q Maxi M spin column with 5 ml diluted Washing Buffer and spin for 5 minutes at 500 × g. We recommend the use of a swing-out rotor to ensure the uniform flow of your sample solution.

C.) Adenovirus purification

15. Load the sample (no more than 20 ml) into the Vivapure Q Maxi M insert and spin the device for 5 min at 500 × g or until the whole volume has passed the membrane. Collect flow-through and repeat the step with residual sample if necessary.

16. Wash the spin column with 18 ml Washing Buffer by spinning for 5 min at 500 × g. Remove flow-through and repeat the wash step once.

17. Using a fresh collection tube, elute adenovirus off the Vivapure Q Maxi M membrane using 1 ml Elution Buffer. Pipette buffer onto the membrane, centrifuge for 30 s at 500 × g and incubate 10 min. Then spin the device for 5 min at 500 × g and collect the adenovirus-containing eluate. Applying a second elution step could improve yield but will dilute virus titre.

Purification protocol – Techniques

D). Optional: Buffer exchange and further concentration

Note:

It is necessary that virus is exchanged into physiological buffer before use in tissue culture or cell based assays, or into generic Storage Buffer for long-term storage at -80°C . Storage Buffers containing glycerine may take considerably longer to concentrate than the original viral eluate solution; prolong centrifuge times and if necessary use cooling at $+4^{\circ}\text{C}$.

18. Transfer eluate to a Vivaspin 20 centrifugal concentrator and add storage|physiological buffer to the concentrate to bring the volume up to 10 ml. Counterbalance the rotor with a second concentrator filled with an equivalent volume of PBS or water. In fixed angle rotors the printed graduations should face away from the centre of the rotor.

19. Centrifuge for 30 min at up to $800 \times g$ in a swing-out rotor, or 25° fixed-angle rotor, with cavities accepting 50 ml conical bottom tubes.

20. Check the volume of viral concentrate remaining in the upper chamber and if necessary centrifuge again and repeat buffer exchange a second time.

⚠ Caution:
Do not reduce the volume to less than 200 μl in order to avoid aggregation and loss of infectivity.

21. Recover the concentrated virus by pipette. Resuspend concentrated virus by gently pipetting up and down a few times before recovery.

22. Determine viral titre. Aliquot accordingly and store virus at -80°C .

Storage buffers for Adenovirus are to be found on page 3 and in the following publication: Hoganson, D. K. et al., Development of a Stable Adenoviral Vector Formulation (2003), Bioprocessing Journal, pp. 43-48.

General information

Typical performance

For a normal yielding vector, 1 × 15 cm culture plates purified using this method should yield a range of up to 1 × 10¹² viral particles (see table 1).

Usage tips

- It is recommended that virus is exchanged into normal physiological buffer before use in tissue culture or cell based assays.

- Aliquot and store virus at -80°C. Once thawed, keep at +4°C and do not re-freeze.

- Virus should remain viable for up to 2 years at -80°C when purified by this procedure.

Table1: Purification results from preparations with certain Ad5 GFP-constructs – depending on individual conditions values may be different.

Purification method	Process time	Eluate	Recovery**	Viral particles
20 ml culture ¹²	1 hour	1 ml	65–70%	1 × 10 ¹¹⁻¹²
500 ml CsCl	12–48 hours	1–2* ml	60–70%	1 × 10 ¹¹⁻¹²

* after dialysis

** before buffer exchange

Ordering information

Ordering Information			Pack Size
VS-AVPQ020	Vivapure AdenoPACK™ 20, 20 ml culture volume		6
VS-AVPQ022	Vivapure AdenoPACK™ 20 RT, 20 ml culture volume*		6
VS-AVPQ101	Vivapure AdenoPACK™ 100, 200 ml culture volume		1
VS-AVPQ102	Vivapure AdenoPACK™ 100 RT, 200 ml culture volume*		1
VS-AVPQ501	Vivapure AdenoPACK™ 500, 500 ml culture volume		1
VS-AVPQ502	Vivapure AdenoPACK™ 500RT, 500 ml culture volume*		1

Additional products in this Kit

VS2041	Vivaspin 20, 100,000 MWCO PES	6
--------	-------------------------------	---

Related Products

5441307H0-00	Sartopore 2 150 0.45–0.2 µm PES	5
17829-K	Minisart Plus 0.45 µm CA + GF	50

* Kit does not contain Benzonase®

Sartorius Lab Instruments GmbH & Co. KG
Otto-Brenner-Strasse 20
37079 Goettingen, Germany

Phone: +49.551.308.0
Fax: +49.551.308.3289
www.sartorius.com

The information and figures contained in these instructions correspond to the version date specified below.

Sartorius reserves the right to make changes to the technology, features, specifications and design of the equipment without notice.

Masculine or feminine forms are used to facilitate legibility in these instructions and always simultaneously denote the other gender as well.

Copyright notice:

This instruction manual, including all of its components, is protected by copyright.

Any use beyond the limits of the copyright law is not permitted without our approval.

This applies in particular to reprinting, translation and editing irrespective of the type of media used.

© Sartorius Germany

Last updated:

11 | 2016