

MS Mini Incubator

MO-MINI

Major Science's Mini Incubator is designed for personal use and small laboratories, saving much of space. The unit features a broad temperature range to meet a variety of microbiology or hematology applications. It is compact and economical, featuring unique features. The housing is all corrosion resistant metal. The glass window of the door offers full visibility to the interior. One stainless steel shelf is included, and can be adjusted to two different levels - additional shelves can be purchased to increase storage capacity. The colorful touch panel presents vivid and intuitive user interface and the backlit LED ensures enough brightness even in dark environment.

The heating structure incorporated in MS Mini Incubator performs evenly heating and precise temperature control. The microprocessor ensures temperature stability within the chamber. A thermometer port on the top of the unit allows the users to monitor temperatures.

In addition, a designed spare port is on the rear of the chamber which allows an electrical cord to be routed out if needed. Optional function control software provides real-time data recording through a PC.



Specification

Cat. No.	MO-MINI
Display	LCD
Temperature Range	Ambient +5° to 70°C
Temperature Accuracy	± 0.2°C @ 37°C in center point
Capacity	17L
Exterior Dimension (WxLxH)	Approx. 12.2"x12"x15" (310x306x380mm)
Interior Dimension (WxLxH)	Approx. 10.3"x10"x10" (261x255x255mm)
Weight	Approx. 28.7 (13kg)
Power	100-240V~, 50/60Hz, 2A
Material	Corrosion resistant metal
PC Connection	USB

Ordering Information


Cat. No	Product Description
MO-MINI	MS Mini Incubator with one stainless steel shelf

Accessories

Cat. No	Product Description
MO-SH2524	250x240mm stainless steel shelf for MS Mini Incubator
MO-MINI-SW	Functional control software for MS Mini Incubator, including a USB cable

Features

- Ideal for microbiology or hematology applications
- Corrosion resistant metal housing
- Door with large viewing area
- Backlit colored touched panel
- One stainless steel shelf is included
- Function control software available

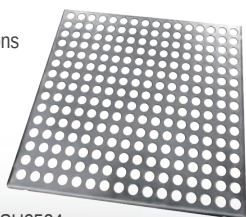

MO-SH2524



* All images are for reference only, actual products might differ from the pictures above.
* Technical specifications subject to change without notice.